



IMPORTANT SAFETY INSTRUCTIONS: All safety and operating instructions should be carefully read before use.

CAUTION: to reduce the risk of electric shock, do not remove cover or open the device. There are no serviceable parts inside.


1. Attach the magnetic charge USB cable to the charging port at the bottom of the massager.

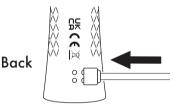


2. Plug the USB end of the cord into the USB port of your computer or other power source. A white light will be blinking on the middle button of the massager indicating battery charging.

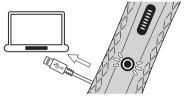

3. The light will become constant when the battery is fully charged. Unplug the USB charging cord. It takes appx 1.5 hour for full charge (depends on the power supply).

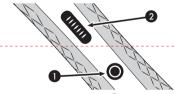

1.Press and hold 'power button ① ' for 2 seconds to turn on the device, LED stays solid=turn it on=go to standby mode. Press again the button ① to Power off.

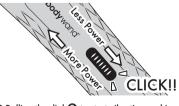

2.Rolling the dial **2** to start vibration and increase or decrease the vibration power.

3. When rolling the dial to the end of less power, the dial will sound 'click!' to make sure turn on or off the vibration. WARNING: For adult use only. For external use only. Remove all scarves, neckties, necklaces, and jewelry before using this massager.

IMPORTANT: Any individual who may be pregnant, has a pacemaker, suffers from diabetes, phlebitis and/or thrombosis, is at an increase risk of developing blood clots, or who has pins/screws/artifical joints or other medical devices implanted in his/her body should consult with a physician before using a massaging device designed for home use.

- Do not reach for an appliance that had fallen into water. Unplug immediately.
 Do not use while bathing or in the shower.
- Do not place or store appliance where it can fall or be pulled into a tub or sink.

 Do not place in or drop into water or other
- 5. Always unplug this appliance from the electrical outlet immediately after use or before cleaning

CARE AND MAINTENANCE

Clean the massager with a soft, damp cloth. DO NOT clean with an alcohol or chemical-based cleaner.

Store at 20 ° - 25 ° Store at cool dry place.

SPECIFICATIONS:

- Power 3.7V/ 2000MAH USB 2.0
- Battery Life: Up to 1 hour use depending on speed and massage platform choosen.
- IPX6 Splash proof for clearing the wand only, can not play under water.

DISCLAIMER:

BODYWAND SHALL NOT BE LIABLE FOR INCIDENTAL OR CONSEQUEN-TIAL DAMAGES RESULTING FROM THE MISUSE OF THIS PRODUCT.



XGEN Products™ Horsham, PA 19044 All Rights Reserved

UK Responsible Person:
ArangaRanga UK LTD, 20-22 Wenlock Road, London, N17 GU, UK.
Phome: +44 2037692960, E-Mail: UKresponsibleperson@arangaranga
EU Responsible Person:
ArangaRanga SL, Paseo de la Castellana 194, Madrid, 28046, Spain.
Phone: +34 910054258, E-Mail: EUresponsibleperson@arangaranga.







